


The logo for SelfEx, consisting of a stylized icon made of three horizontal bars in blue, purple, and red, followed by the text "SelfEx" in a large, bold, sans-serif font.

Peripheral Self-Expandable Stent System

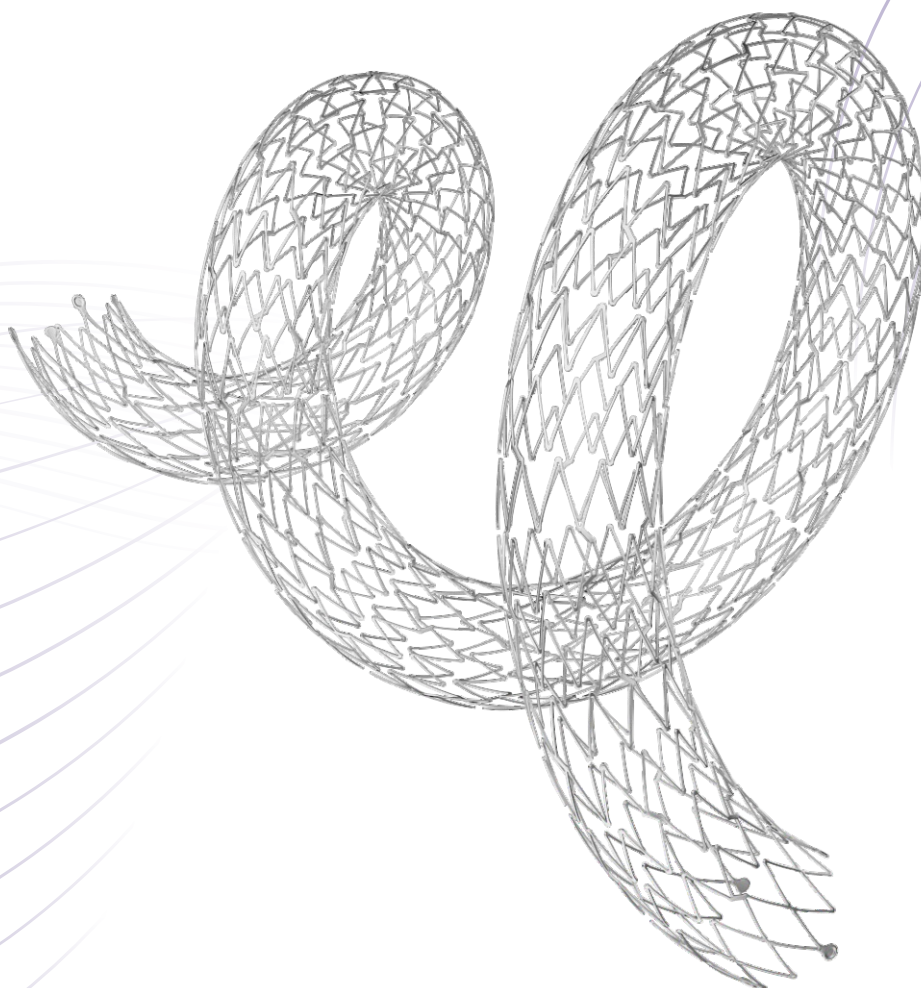

ENGINEERED FOR EXCELLENCE

TRUSTED PERFORMANCE

- ≡ Stent Architecture with nitinol material provides shape memory allowing the stent to be easily deployed
- ≡ Multifunctional braided catheter for optimal balance between pushability, trackability and flexibility

SUPERIOR FLEXIBILITY



- ≡ Hybrid stent design provides multidirectional flexibility
- ≡ Crafted with closed cell ends guarantees deployment stability
- ≡ Open cell centered provides unparalleled flexibility & heightened fracture resistance

RADIOPACITY

- ≡ Engineered for enhanced visibility under fluoroscopy
- ≡ Attached tantalum marker increases visibility and precise deployment

RADIAL STRENGTH

- ≡ Ability to resist compression & maintain luminal gain
- ≡ Dimensional design changes facilitated harmonisation of radial force across the range of product diameters
- ≡ SelfEx's[®] shape memory properties enable the stent to apply a gradual and persistent chronic outward forces, promoting healing without causing excessive stress or damage to the vessel wall

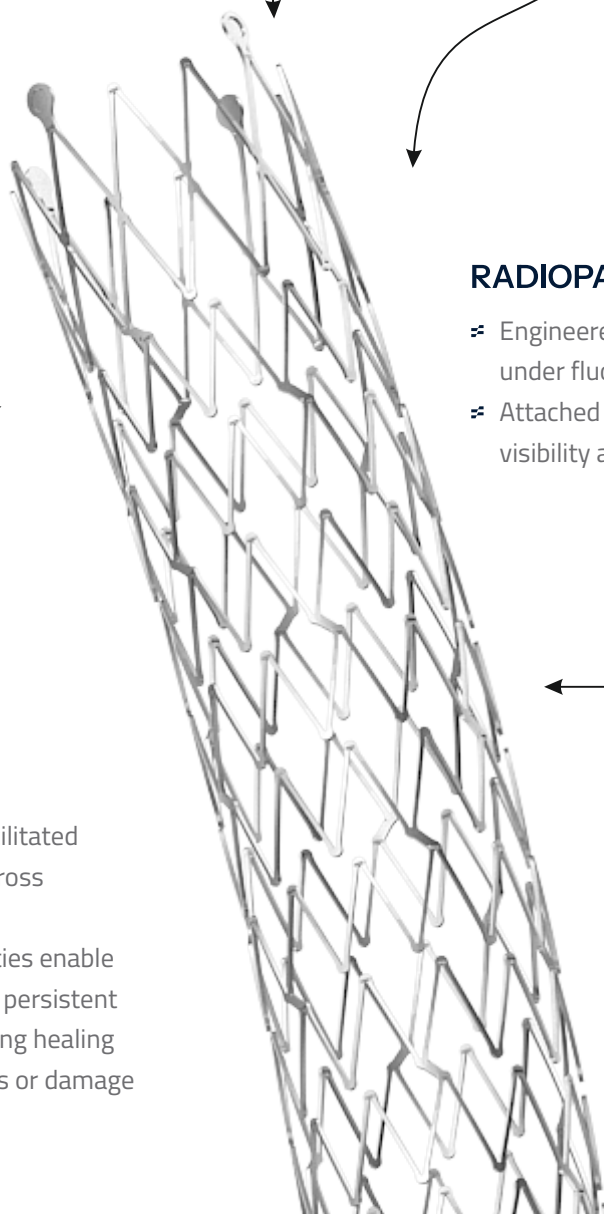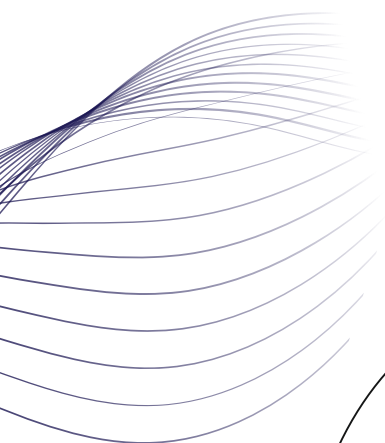

EXCELLENT DEPLOYMENT

- ✦ One handed stent release with ergonomically designed handle
- ✦ Offers ease of use with accurate deployment



STENT

Stent Material	Nickel Titanium Alloy
Design Features	Hybrid Cell stent design
Stent Diameters (mm)	4.00, 5.00, 6.00, 7.00, 8.00, 9.00, 10.00
Stent Lengths (mm)	20, 30, 40, 60, 80, 100, 120, 150, 180, 200, 220

DELIVERY SYSTEM

Delivery System	Over-the-wire (OTW)
Shaft Outer Diameter	6F
Shaft Marker Bands	Platinum – Iridium bands
Usable Catheter Length	120 cm
Recommended Minimum Sheath / Introducer	6F
Guide wire Compatibility	0.035" (0.89 mm)

ORDERING INFORMATION

Nominal Stent Length (mm)	Unconstrained Stent Diameter						
	4.0 mm	5.0 mm	6.0 mm	7.0 mm	8.0 mm	9.0 mm	10.0 mm
20	SE0400020	SE0500020	SE0600020	SE0700020	SE0800020	SE0900020	SE1000020
30	Se0400030	SE0500030	SE0600030	SE0700030	SE0800030	SE0900030	SE1000030
40	SE0400040	SE0500040	SE0600040	SE0700040	SE0800040	SE0900040	SE1000040
60	Se0400060	SE0500060	SE0600060	SE0700060	SE0800060	SE0900060	SE1000060
80	SE0400080	SE0500080	SE0600080	SE0700080	SE0800080	SE0900080	SE1000080
100	Se0400100	SE0500100	SE0600100	SE0700100	SE0800100	SE0900100	SE1000100
120	SE0400120	SE0500120	SE0600120	SE0700120	SE0800120	SE0900120	SE1000120
150	SE0400150	SE0500150	SE0600150	SE0700150	SE0800150	SE0900150	SE1000150
180	SE0400180	SE0500180	SE0600180	SE0700180	SE0800180	SE0900180	SE1000180
200	SE0400200	SE0500200	SE0600200	SE0700200	SE0800200	SE0900200	SE1000200
220	SE0400220	SE0500220	SE0600220	SE0700220	SE0800220	SE0900220	SE1000220

Corporate Office - Sahajanand Laser Technology Ltd.

E-30, Electronic Estate, G.I.D.C., Sector-26, Gandhinagar - 382 028, Gujarat, India.

Tel : +91-79-6170-1616 | +91-90-9904-0080

Email : mkt@SLTL.com

www.SLTLmedical.com

Manufacturing Site - Sahajanand Laser Technology Ltd.

A-8, Electronic Estate, G.I.D.C., Sector - 25, Gandhinagar - 382 016, Gujarat, India.

Tel : +91-90-9904-0090

SE/CAT/00



Caution: Indication, Contraindication, Warnings, Precautions and other instructions for use can be found in the product labeling supplied with each device.

SelfEx® is registered trademark of Sahajanand Laser Technology Ltd. © 2024 SLTL Group. All right reserved.

SLTM14PC10SE24DT