


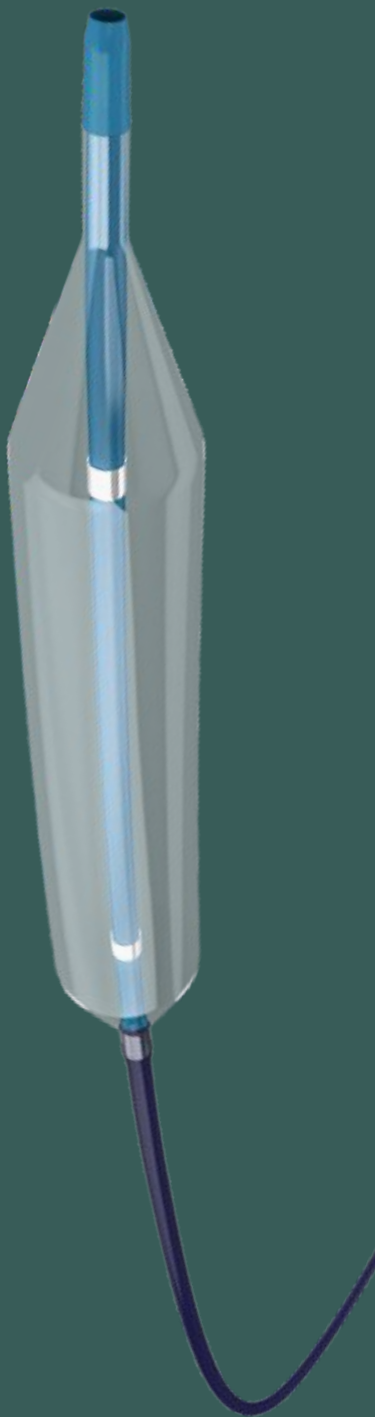


Vector[®] PTA

PTA Balloon Catheter

**BUILD FOR LONGER
STRONGER DILATION**

Vector[®] PTA

- Engineered to enhance performance at every stage of your procedure, combining precision and reliability.
- Designed to keep you ahead in every case with consistent, high-level performance.

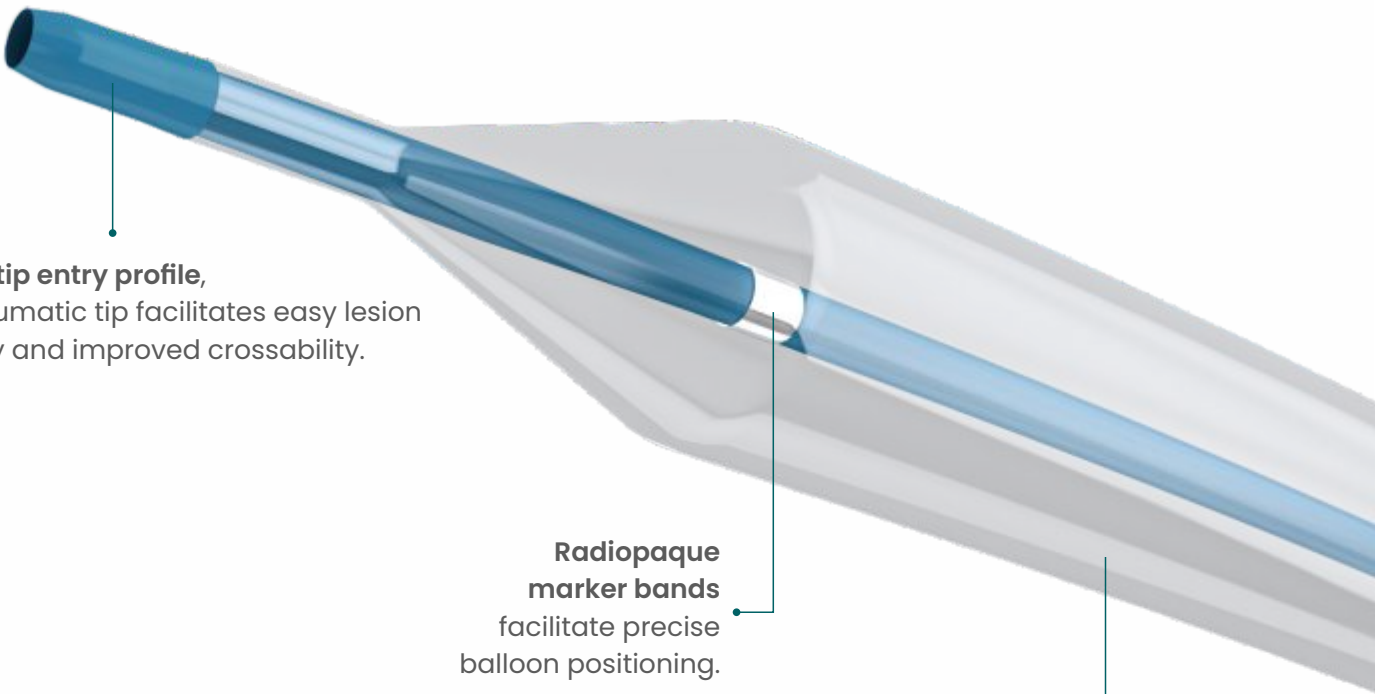


Low tip entry profile,
Atraumatic tip facilitates easy lesion
entry and improved crossability.

**Radiopaque
marker bands**
facilitate precise
balloon positioning.


Six-Fold balloon
enables smooth deflation and
easy catheter withdrawal.

Broad Balloon Diameter Range

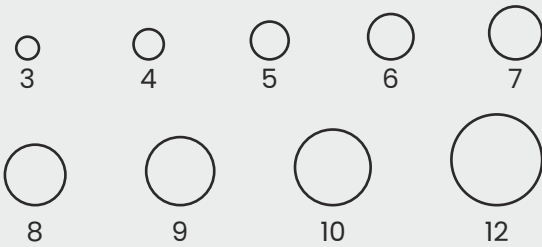

Coaxial Catheter Shaft

Supports rapid inflation, deflation, trackability, and pushability.

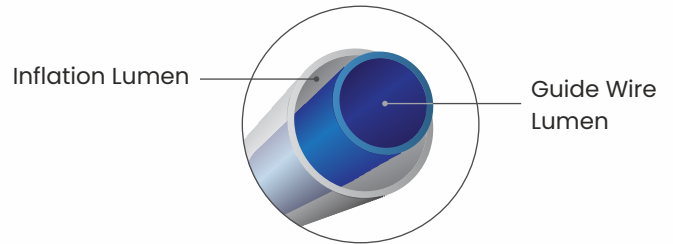

TECHNICAL SPECIFICATIONS

Balloon Diameter (mm)	3, 4, 5, 6, 7, 8, 9, 10, 12
Balloon Length (mm)	20, 30, 40, 60, 80, 100, 120, 150, 180, 200
Radiopaque Marker	2-Platinum Iridium Marker
Effective Catheter Length	135mm
Guide Wire Compatibility	0.035"
Rated Burst Pressure(RBP)	14-18 atm
Nominal Pressure (NP)	8-10 atm

Nylon balloon material

enables great trackability and high pressure resistance to facilitate crossing of tight lesion.

Hydrophilic-coated shaft

ensures excellent trackability and smooth navigation through complex and tortuous vessels.

INDICATIONS

PERIPHERAL ARTERY DISEASE

- Iliac Artery
- Popliteal Artery
- Femoral Artery
- Renal Artery ...

COMPLIANCE CHART

Diameter (mm)	Nominal Pressure (atm)	Rated Burst Pressure (atm)	Inroducer Sheath (Fr)	Guide Wire Compatibility
3	10	18	5	0.035"
4	10	18	5	
5	10	18	5	
6	10	18	5	
7	10	18	6	
8	10	18	6	
9	8	18	6	
10	8	14	6	
12	10	14	8	

ORDERING INFORMATION

Balloon Diameter (mm)	Balloon Length (mm)									
	20	30	40	60	80	100	120	150	180	200
3	VPT35A03020	VPT35A03030	VPT35A03040	VPT35A03060	VPT35A03080	VPT35A03100	VPT35A03120	VPT35A03150	VPT35A03180	VPT35A03200
4	VPT35A04020	VPT35A04030	VPT35A04040	VPT35A04060	VPT35A04080	VPT35A04100	VPT35A04120	VPT35A04150	VPT35A04180	VPT35A04200
5	VPT35A05020	VPT35A05030	VPT35A05040	VPT35A05060	VPT35A05080	VPT35A05100	VPT35A05120	VPT35A05150	VPT35A05180	VPT35A05200
6	VPT35A06020	VPT35A06030	VPT35A06040	VPT35A06060	VPT35A06080	VPT35A06100	VPT35A06120	VPT35A06150	VPT35A06180	VPT35A06200
7	VPT35A07020	VPT35A07030	VPT35A07040	VPT35A07060	VPT35A07080	VPT35A07100	VPT35A07120	VPT35A07150	VPT35A07180	VPT35A07200
8	VPT35A08020	VPT35A08030	VPT35A08040	VPT35A08060	VPT35A08080	VPT35A08100	VPT35A08120	VPT35A08150	VPT35A08180	VPT35A08200
9	VPT35A09020	VPT35A09030	VPT35A09040	VPT35A09060						
10	VPT35A10020	VPT35A10030	VPT35A10040	VPT35A10060						
12	VPT35A12020	VPT35A12030	VPT35A12040	VPT35A12060						

Corporate Office - Sahajanand Laser Technology Ltd.

E-30, Electronic Estate, G.I.D.C., Sector-26,
Gandhinagar - 382 028, Gujarat, India.
Tel : +91-79-6170-1616
Cell : +91-99-2503-6495
Email : mkt@SLTL.com

Manufacturing Site - Sahajanand Laser Technology Ltd.

A-8, Electronic Estate, G.I.D.C., Sector - 25,
Gandhinagar - 382 016, Gujarat, India.
Cell : +91-99-2503-6495

www.SLTLmedical.com
www.SLTL.com



CAUTION: Indication, Contraindication, Warnings, Precautions and other instructions for use can be found in the product labeling supplied with each device.