


zephyr[®]xxl

Cobalt Chromium Peripheral Stent



A Reliable Solution for Extra-Large Vessels

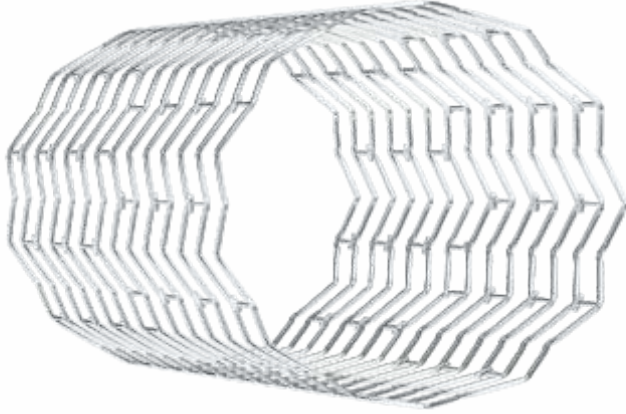

CoCr Stent for Extra Large Vessels:

- 12 crown stent design with 'C' & 'S' shaped polylinks
- Precise positioning and stability in complex cases.

Extra-Large Size & Coverage:

- Meticulously crafted to support vessel upto 30mm.
- High Radial Force maintains vessel patency and optimal lumen diameter.

Atraumatic Edges & Flexibility:

- Adapts to vessel contours for smooth deployment and optimal fit, reducing complications.
- Atraumatic edges with C&S connectors for uniform expansion and consistent outcomes.
- Harmonic strut expansion ensures flexible vessel integration.

Clinical Indications:

- Coarctation of Aorta
- Right ventricular Outflow Tract (RVOT)
- Conduit Stenosis

Technical Specifications:

Stent Material Cobalt Chromium	Stent Design Zigzag Hybrid	Unexpanded Stent Diameter Ø 4.5 mm	Minimum Crimping Profile 10F + 1F
Stent Recoil ◀ 3%	Configuration 12 Crowns	Expanded Stent Diameter Ø 26mm- 30mm	Over Expansion Range Up to 32 mm
Foreshortening ◀ 5%	Strut Thickness 0.009"	Recommended Sheath of Balloon + 2F of Stent	Guide Wire Compatibility 0.035"

Ordering Information:

	Length (mm)										
	24	36	42	50	54	60	66	72	80	85	100
Product Codes	ZPXXL24	ZPXXL36	ZPXXL42	ZPXXL50	ZPXXL54	ZPXXL60	ZPXXL66	ZPXXL72	ZPXXL80	ZPXXL85	ZPXXL100

Corporate Office - Sahajanand Laser Technology Ltd.

E-30, Electronic Estate, G.I.D.C., Sector-26,
Gandhinagar - 382 028, Gujarat, India.
Tel : +91-79-6170-1616
Cell : +91-99-250-36495
Email : mkt@SLTL.com

Manufacturing Site - Sahajanand Laser Technology Ltd.

A-8, Electronic Estate, G.I.D.C., Sector - 25,
Gandhinagar - 382 016, Gujarat, India.
Cell : +91-99-250-36495

www.SLTLmedical.com
www.SLTL.com