


Understanding the Art of Stenting

EasyRap

Stainless Steel Rapamycin (Sirolimus)
Eluting Coronary Stent System

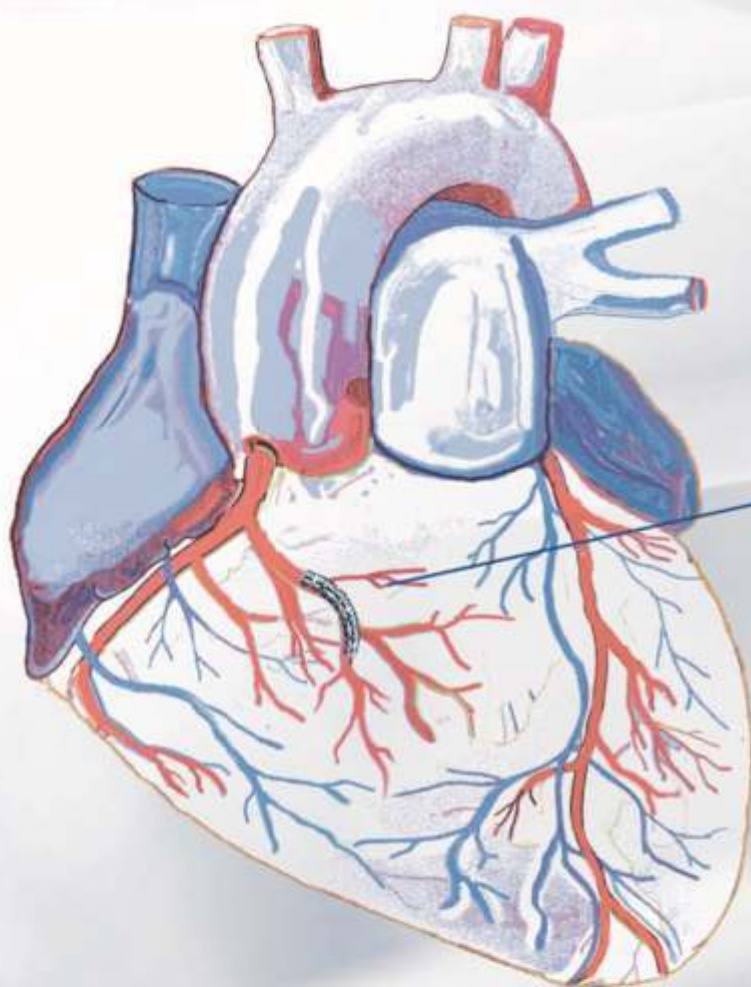

based on the
clinically proven

Easyflex
stent design

EasyRap

Stainless Steel Rapamycin (Sirolimus) Eluting Coronary Stent System


EasyRap

Stainless Steel Rapamycin (Sirolimus) Eluting Coronary Stent System


Biosteel coating technology
Material: Stainless steel
DLC (diamond like carbon coated) covered with polymer & drug
Therapeutic release within 28 days
High mechanical stability of the coating

Technical Specifications

Polymer	Inert polymer compound / fully resorbable
Release kinetic	Slow release within 10-28 days
Stent surface finishing	DLC PII (plasma ion implantation)

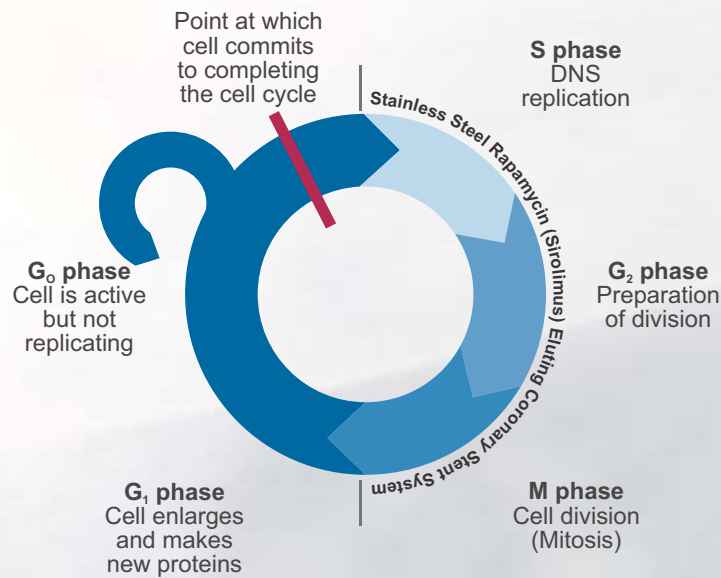

Two opposite articulations between every stent segment guarantee multidirectional flexibility

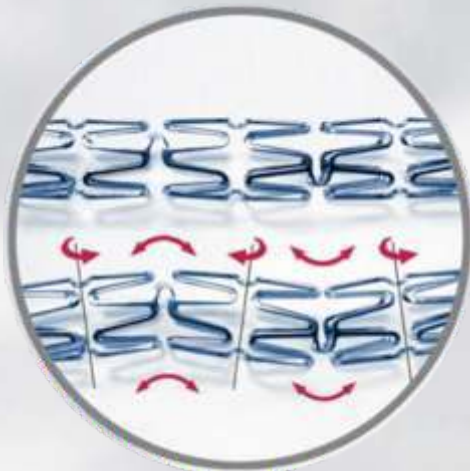

Drug
Rapamycin

Target Dose
1µg/mm²

Polymer Type
Biodegradable

Strut Thickness
90µ

Coating thickness
5-6µ

Balloon Overhang
<0.5mm

Usable Catheter Length
145cm

Guide wire Compatibility
0.014"

Guide catheter Compatibility
5F

MRI Compatibility
MR Conditional

Mean Recoil
2.3%

Mean Foreshortening
0.6%

Crossing Profile
0.97 mm (3mm)

Therapeutic Release
Slow release up to 28 days

Stent material
316LVM Stainless Steel

EasyRap

Stainless Steel Rapamycin (Sirolimus) Eluting Coronary Stent System

CATALOG NO :

Length	Diameter (Ø)							
	2.00mm	2.25mm	2.50mm	2.75mm	3.00mm	3.50mm	4.00mm	4.50mm
07mm	ER 0720	ER 0722	ER 0725	ER 0727	ER 0730	ER 0735	ER 0740	ER 0745
10mm	ER 1020	ER 1022	ER 1025	ER 1027	ER 1030	ER 1035	ER 1040	ER 1045
13mm	ER 1320	ER 1322	ER 1325	ER 1327	ER 1330	ER 1335	ER 1340	ER 1345
15mm	ER 1520	ER 1522	ER 1525	ER 1527	ER 1530	ER 1535	ER 1540	ER 1545
17mm	ER 1720	ER 1722	ER 1725	ER 1727	ER 1730	ER 1735	ER 1740	ER 1745
20mm	ER 2020	ER 2022	ER 2025	ER 2027	ER 2030	ER 2035	ER 2040	ER 2045
24mm	ER 2420	ER 2422	ER 2425	ER 2427	ER 2430	ER 2435	ER 2440	ER 2445
28mm	ER 2820	ER 2822	ER 2825	ER 2827	ER 2830	ER 2835	ER 2840	ER 2845
33mm	ER 3320	ER 3322	ER 3325	ER 3327	ER 3330	ER 3335	ER 3340	ER 3345
38mm	ER 3820	ER 3822	ER 3825	ER 3827	ER 3830	ER 3835	ER 3840	ER 3845
42mm	ER 4220	ER 4222	ER 4225	ER 4227	ER 4230	ER 4235	ER 4240	ER 4245
45mm	ER 4520	ER 4522	ER 4525	ER 4527	ER 4530	ER 4535	ER 4540	ER 4545

International Head Quarter - Sahajanand Laser Technology Ltd.

E-30, G.I.D.C., Electronic Estate, Sector-26,

Gandhinagar - 382 028. (Gujarat) India.

Tel. : +91 - 79 - 2324 1379 / 80,

Fax : +91 - 79 - 2324 1465

Email : mkt@SLTL.com

Log on to :

www.SLTLMedical.com

www.SLTL.com

India Marketing - Sahajanand Laser Technology Ltd.

41, Newyork Tower-A, Nr. Thaltej Cross Road,

S.G. Highway, Ahmedabad - 380 054. Gujarat, India.

Tel.: +91-79-2685 4537 / 38 , Fax: +91-79-2685 4539

Manufacturing - Sahajanand Laser Technology Ltd.

A-8, G.I.D.C., Electronic Estate, Sector-25,

Gandhinagar - 382 016. (Gujarat) India.

Tel.: +91 - 79 - 2328 7462, Fax: +91 - 79 - 2324 1465



CAUTION: Indication, Contraindication, Warnings, Precautions and other instructions for use can be found in the product labeling supplied with each device.
EasyRap is a registered trademark of Sahajanand Laser Technology Ltd. © 2020 SLTL. All right reserved.